



CHERISH. EVERY. BREATH.

iGo<sup>®</sup>2

### Drive DeVilbiss iGo<sup>®</sup>2 Portable Oxygen Concentrator (125D)

- Automatically adjusts oxygen settings – so you don't have to!
- Weighs less than 5 lbs.!
- One of the quietest units on the market
- Easy to carry in its convenient carry bag
- Exceptionally durable construction for dependable use


SmartDose<sup>®</sup>  
TECHNOLOGY



Monitors breathing patterns, and adjusts dosing to match activity at every breath.

CELEBRATING  
**130**  
YEARS

For more than 130 years, DeVilbiss products have paved the way for the respiratory industry. Here's another product they'll love from the name you can trust.

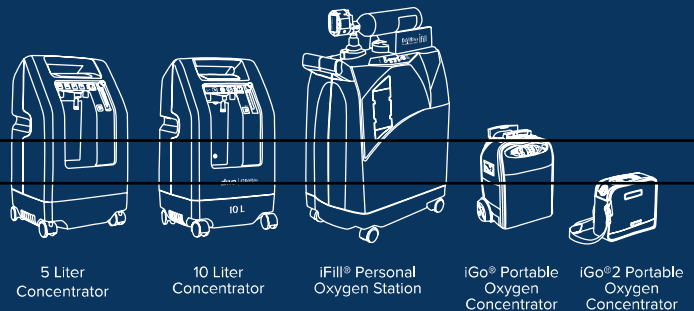

5 Liter Concentrator

10 Liter Concentrator

iFill<sup>®</sup> Personal Oxygen Station

iGo<sup>®</sup> Portable Oxygen Concentrator

iGo<sup>®</sup>2 Portable Oxygen Concentrator

**WEIGHS LESS THAN 5 LBS.!**

## Why iGo2 POC?

The first-of-its-kind iGo2 portable oxygen concentrator uses patented SmartDose® Technology to automatically adjust your oxygen setting.

**How does the SmartDose® Technology work?** The technology detects your breath rate and will deliver more oxygen when it increases, ensuring that your oxygen level stays consistent.

You'll no longer have to manually adjust oxygen settings and will gain the freedom and confidence to live a more active lifestyle!

### Features:

- Weighs less than 5 lbs. with battery installed
- Easy-to-read and -access LCD screen
- Readily accessible battery pack
- Fits in a convenient carry bag with pockets for extra batteries
- Choice of delivery modes – SmartDose Technology or standard PulseDose® mode
- More sensitive conserver trigger with the ability to detect your breathing at a lower pressure than other concentrators— you receive the oxygen you need without delay!
- Protective, rugged overmold provides an extra layer of protection against drops and damage
- Federal Aviation Administration (FAA) compliant for air travel

### Accessories



USA Power Cord  
DV51D-606



DC Power Cord  
DV6X-619



Carrying Case  
125D-670



iGo2 Charging Station  
125CH-613  
(Battery not included)



Spare Battery  
125D-613

### Specs

Dimensions	8.4"(H) x 3.5"(W) x 8.6"(L) (21.3cm x 8.9cm x 21.8cm)
Weight	4.95 ± 0.1 lbs.
Oxygen Output	90% + 4/-3 Over Entire Operating Range
Oxygen Delivery Settings	1 to 5
Maximum O <sub>2</sub> Production	1014 mL
Measured Sound Level	41 dBa
Approximate Battery Run Time	3.5 Hours at Setting 2



#### Drive DeVilbiss International

DeVilbiss Healthcare GmbH  
Kamenzer Straße 3  
68309 Mannheim, Germany

UK +44 (0) 121 / 521 31 40  
DE +49 (0) 621 / 1 78 98 0  
kontakt@devilbisshc.com



www.drivedevilbiss-int.com

