

CPAP | Total Face Mask
Type 837



Total Face Mask
(For Adult)
#13843 Small - S
#13844 Medium - M
#13845 Large - L



Total Face Mask
(For Pediatric)
#13831 X Small
#13832 Small



Total Face Mask
(For Pediatric)
#13830 XX Small



 **HSINER**[®]
Touching lives, Improving life.

CPAP Total Face Mask



 **HSINER**[®]
Touching lives, Improving life.



HSINER Co., Ltd.

No. 312, Jhongshan Rd., Shengang Dist., Taichung City 429, Taiwan.


Tel: 886-4-2515-2480 Fax: 886-4-2515-2482

<http://www.hsiner.com> E-mail: sales@hsiner.com




CPAP Total Face Mask

Type 837

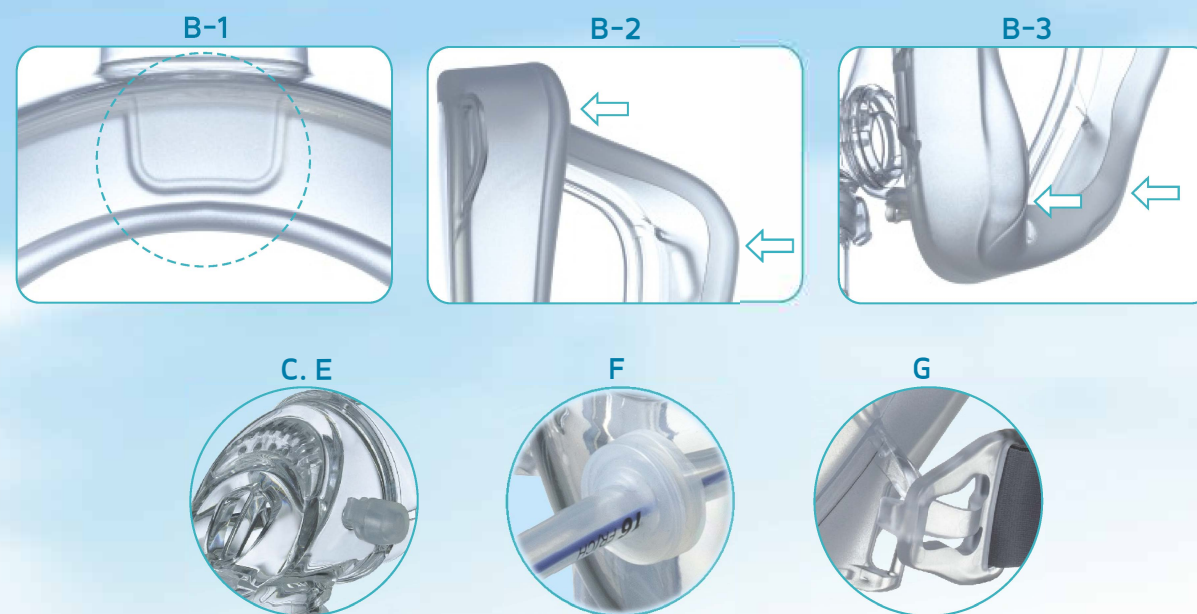

- A. The transparent cover creates a crystal clear view for the user.
- B. Silicone total full face cushion fits just right and firmly on the whole face, has less pressure for the user.
 - B-1 : Strengthened stiffness of the cushion decreases leakage, without tightening the headgear, easing the pressure to the user.
 - B-2 : Flexible mask frame fits different facial contours, providing better sealing and comfort.
 - B-3 : Ergonomic mask curve allows user to open mouth with the mask on.
- C. Anti-asphyxia valve, automatically open when the generator pressure fails, enabling patients to breathe by themselves.
- D. Quick release adaptor, press make user release / attach the tubing easily during temporary leaving.
- E. Sampling port allows measuring pressure or ETCO2.
- F. With nasogastric (NG) port.
- G. Quick release buckles, easy for user to release and take off the mask after sleep.
- H. 360° rotational elbow providing a comfortable tube position.
- I. Applicable to the users who have facial hair and general facial irregularity.

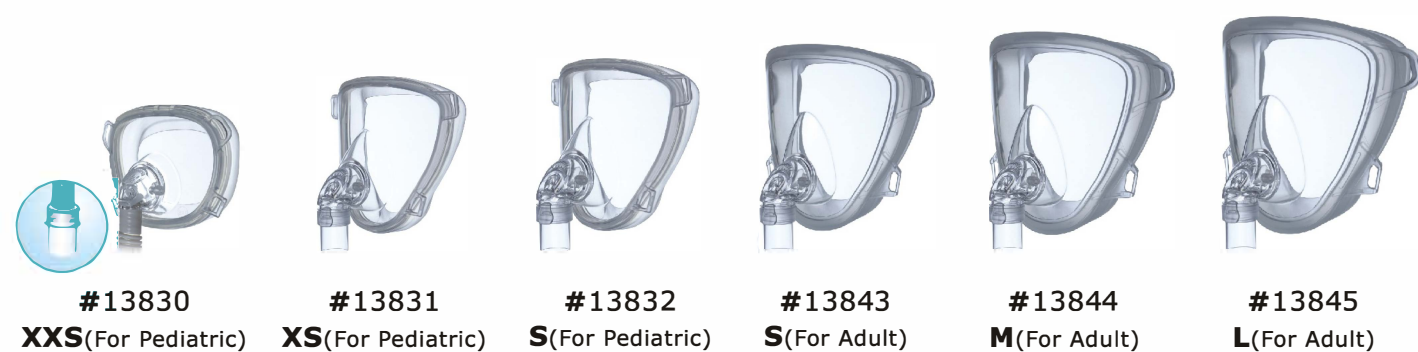

Optional Elbow Connectors

