


F&P Nivairo+™ with ErgoFit™ 2 Headgear
Take care of the pressure



Nivairo+ helps you address the **NIV points of care**

Noninvasive ventilation (NIV) therapy can be stressful for both the caregiver and patient. There are three interdependent points of care in NIV therapy which typically challenge clinicians daily.

Nivairo+ helps take care of the pressure along the NIV journey. With ErgoFit 2, it is designed to help you achieve a personalized fit for a range of patient head sizes.

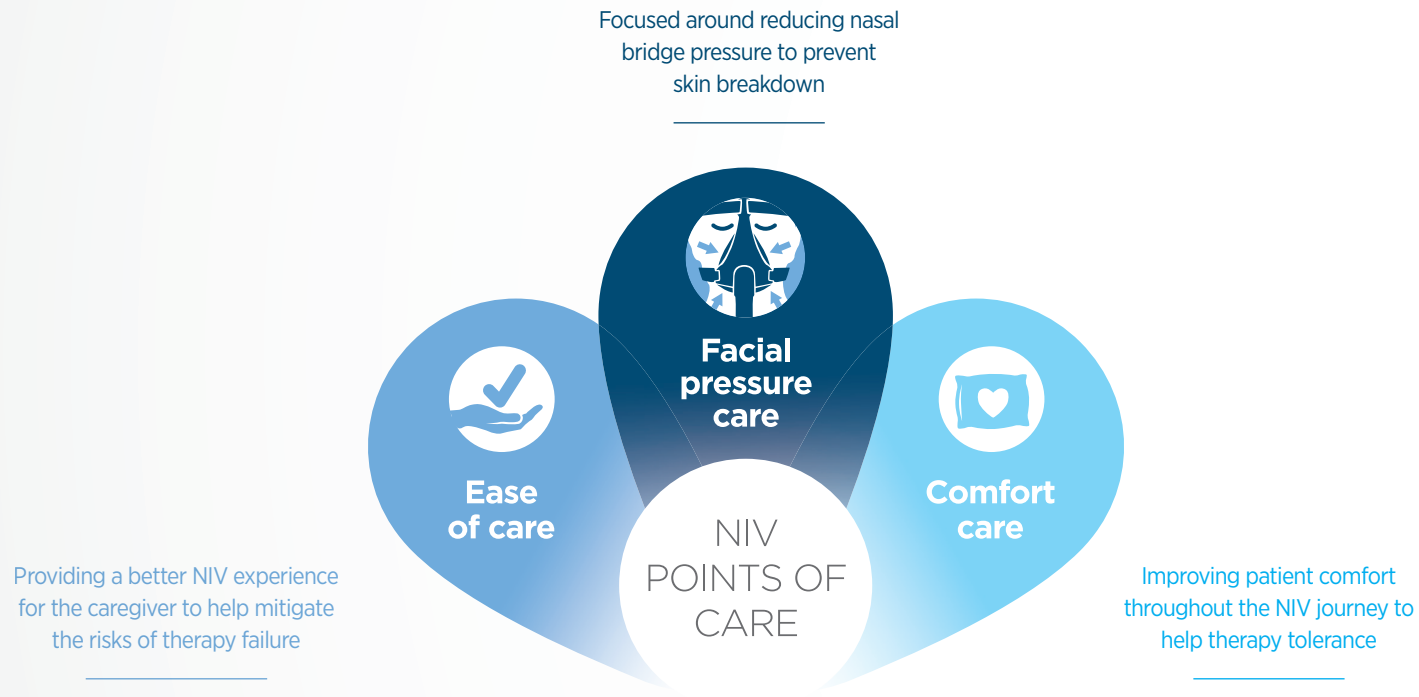



NIV masks have been cited as one of the leading causes of pressure ulcers in the ICU.¹

A hospital-acquired pressure injury can cost a hospital in the US from \$20,900 to \$151,700 per pressure ulcer.²

Nivairo+ with ErgoFit2 is designed to provide a soft seal that helps reduce excessive pressure on the bridge of the nose.



Patient discomfort may be due to irritation from leakage, incorrect mask sizing or over-tightening of the headgear.

Patient discomfort remains responsible for up to 33% of NIV failure.³

Nivairo+ with ErgoFit 2 is designed to help deliver a more personalised fit for your patients to help increase comfort and tolerance.



It can be challenging to fit a mask precisely. There are an infinite number of face shapes and sizes to contend with. This can make fitting a mask challenging.

Nivairo+ with ErgoFit 2 is designed to provide an easy sizing, fitting and readjusting experience to help you provide quality NIV therapy.



Nivairo+ Features

vvvv



RollFit™ Auto-Adjusting Seal

Designed to help fit the unique contours of the nasal bridge



TubeFit™ Tube Sealing Zones

Designed to help maintain an effective seal around nasogastric and orogastric tubes



SizeMe™ Mask Sizing Guide

Simple, colour coded design to help choose the right mask size

ErgoFit 2 Headgear Features



Adjustable Crown Strap

Designed to help you provide a more personalized fit for your patient



Stronger Headgear

Designed to help maintain comfort and structural integrity as you fit and readjust your patients



Two Headgear Sizes

Designed to help accommodate a range of patient head sizes according to their mask size

F&P Nivairo+ suits a variety of NIV patients

Choose from three configurations to suit your unit



RT055

Anti-asphyxiation valve
elbow for single-limb circuits



RT056

Standard elbow
for dual-limb circuits



RT057

Vented with anti-asphyxiation
valve for single-limb circuits



4 Sizes

Four mask sizes available for all mask
configurations including an extra small

F&P Optiflow™ as an NIV rest therapy

F&P Nivairo+ masks are compatible with F&P humidifiers to allow for a seamless NIV experience.

Optiflow is also a rest therapy option to enable patients to take a comfortable break from NIV.⁴

| F&P Nivairo+ Hospital Full Face Masks | Product Order Codes | | |
|---|---|-------------------------------|--|
| | Anti-Asphyxiation Valve Nonvented Version | Standard Elbow | Anti-Asphyxiation Valve Vented Version |
| F&P Nivairo+ X-Small | RT055XS | RT056XS | RT057XS |
| F&P Nivairo+ Small | RT055S | RT056S | RT057S |
| F&P Nivairo+ Medium | RT055M | RT056M | RT057M |
| F&P Nivairo+ Large | RT055L | RT056L | RT057L |
| Mask Type | Nonvented with anti-asphyxiation valve | Nonvented with standard elbow | Vented with anti-asphyxiation valve |
| Operating Pressures | 4 - 25 cmH ₂ O | | |
| Box Quantity | 10 | | |
| Intended Use | Single-patient use; maximum 14 days | | |
| Materials | SEBS Thermoplastic Elastomer, Silicone, Acetal, Polypropylene, Polycarbonate, Polyester Fabric, Polyester Polyurethane Foam, Nylon, Nylon Spandex | | |
| Manufacturing Mode | Non invasive device; medically clean | | |
| Destruction | Incineration, or according to hospital protocol | | |
| Classification | RT055 XS/S/M/L, RT056 XS/S/M/L, RT057XS/S/M/L, USA - Class 2; AU - Class IIa; EU - Class IIa; China - Class 2, Japan - Class I | | |
| Notified Body | TÜV Product Services GmbH | | |
| Circuit Setup | Single-limb and exhalation port | Dual limb | Single limb |
| Recommended F&P 850 System Breathing Circuits | RT319/RT219 | RT380 | RT319/RT219 |
| Recommended F&P 950 System Breathing Circuits | 950A60/950A61 | 950A81/950A82 | 950A60/950A61 |

618122 REV C © 2019 Fisher & Paykel Healthcare Limited

F&P Nivairo, Nivairo+, RollFit, TubeFit, ErgoFit, SizeMe, F&P Optiflow, F&P 850, F&P 950, Intelligent Seal and Sensitive Care are trademarks of Fisher & Paykel Healthcare Limited. For patent information, see www.fphcare.com/ip

Products in this brochure may not be available in all regions.

References. 1. Clay, P. J Ped Nursing 2018;41, 77-79. 2. Agency for Healthcare Research and Quality 2012. 3. Carron, M. Brit J Anaesthesia 2013;110(6), 896-914. 4. Spoletini G. J CritCare 2018; 48, 418-425.